Preva DC

WHERE TECHNOLOGY MEETS SIMPLICITY...PERFECTLY



Preva

In 2002, the Progeny[®] design team set out to create the perfect intraoral x-ray. And they did...first with the reliable JB-70 AC system, followed by the extraordinary, industry-leading, high frequency system, the Preva DC.

Preva is the convergence of the latest DC technology, microprocessor advances and userfocused features. Designed to provide flawless imaging and enhanced workflow, Preva is the perfect foundation for every practice—film-based, PSP and digital—and offers your practice reliable performance today, tomorrow and for many years to come.

Now there are two choices in imaging excellence: Preva the perfect foundation for every practice and every receptor; and Preva Plus, the DC intraoral x-ray with its own integrated digital sensor. Only Preva into one

progeny

60kv 7mA

0100s

Preva

V

Ŷ.

A) 0

 (\triangleright)

BAA

packs so many powerful features attractive, reliable DC system.

The smallest focal spot equals the sharpest images


Preva's 0.4mm focal spot produces visibly higher resolution images whether you're using film or digital sensors. Details are sharper and more defined. Small but important differences are easier to spot. Pathologies are easier to identify.



Traditional 0.7mm focal spot



Preva 0.4mm focal spot

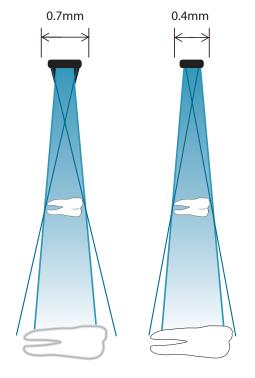

Adjustable techniques provide complete control

With Preva, you decide how you want to see your images. A choice of selectable options allows you to determine the techniques you prefer to control gray scale, density and contrast according to clinical or procedural needs. And, only Preva combines a 0.4mm focal spot with such a wide range of kV, mA and time choices for consistently superior diagnostic images.

- n kV 60, 65 and 70
- n mA 4, 5, 6 and 7
- Time 0.02 seconds to 2 seconds



All the benefits of Preva, PLUS a high resolution digital sensor always handy, always ready.

progent

Sensor stores safely and securely on arm

Integrated efficiency

Preva Plus takes digital imaging to the next level of convenience with its innovatively integrated ClearVision[®] sensor. All cabling and electronics are self-contained and hidden away inside the Preva Plus, helping to keep your operatory well-organized and cable-free.

Sensors that are easily misplaced or damaged can now be safely stored within easy reach. A simple USB connection and a specially designed receptacle right on the arm makes it easy to be ready for your next patient.

> The latest USB technology provides "plug-and-play" installation.

All the details, no wait

Progeny's[®] ClearVision[®] sensor produces crisply detailed digital images ready for viewing on your monitor within seconds. Images can be captured and viewed with Progeny[®] Imaging, and most popular Practice Management and Imaging Programs.

For a complete list of software compatibility, contact Midmark.



Enhanced wo



Stable, drift-free positioning

Preva is known for its exceptional arm mechanics. Feel the exceptional balance and effortless movement. Preva's proprietary braking system contains both an outboard and an inboard brake to assure that the head stays exactly where you position it without sagging or drifting. Preva has more adjustment points than any other system and a distinctive handle design for fast onehanded positioning.





Only Preva offers a dual arm braking system

orkflow

Tailor-made for optimum efficiency

Preva's exclusive LCD display with its large, bright screen displays technique information clearly from any angle. Easy to understand icons quickly identify tooth, receptor type and patient size. 30 pre-programmed technique settings streamline setup and preparation for every procedure and every receptor.

> 56", 66", 76" and 82"* arm lengths available *Requires mounting plate



Optimize your Preva

Only Preva delivers menu options that allows you to instantly optimize your x-ray. Use the built-in menu accessed from the control panel to change preset techniques or select a specific sensor. Progeny's[®] unique Optimization Menu puts many options within your reach—simply.



progeny

Preva

B

Ŷ


θ

mΑ

M

Ð

Progeny[®]—the leader in intraoral x-ray systems in the U.S. and the preferred choice of dentists around the world.

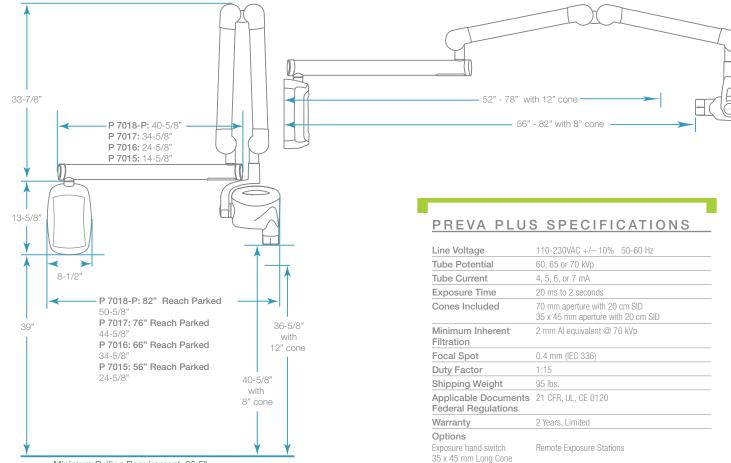

Minimum Ceiling Requirement: 86.5"

Cabinet Pass Through 66" Arm: 31" w x 47.5" d 76" Arm: 41" w x 47.5" d

Cabinet Pass Through (Alternate) 66" Arm: 26.75" w × 47.5" d 76" Arm: 36.75" w × 47.5" d

Mobile Base: 26.9" w x 36" | Height: 76.9" Arm reach: 44"

midmark.com

Midmark Corp. 675 Heathrow Drive Lincolnshire, II 60069 847-415-9800 Toll-free 888-924-3800

©Midmark 2013 All Rights Reserved Printed in the USA ML-00041 Rev. D

Computer Requirements

Choice of 56", 66", 76", or 82" (total reach) and Mobile

Size 1

90 cm

Two stud mounting plate

Film Size equivalent

Computer Interface

Sensor

Dimensions

Active Area

Sensor Cable

Please see www.progenydental.com for operating system and hardware requirements

Metal stud mounting plate

37 x 24 x 5.5 mm

High Speed USB 2.0

30 x 20 mm

Contact your Midmark representative for additional system requirement information.



Size 2

90 cm

43 x 32 x 5.5 mm

High Speed USB 2.0

36 x 26 mm