

ULTRA CLEAN OILLESS AIR COMPRESSOR



- **COMPRESSOR DESIGN OFFERS THE QUIETEST COMPRESSORS IN THE INDUSTRY**
- **ALL AIR RECEIVER TANKS HAVE INTERNAL EPOXY COATING APPLIED**
- **THE AIR IS DRYED BEFORE THE TANK TO PREVENT BACTERIA BUILD UP IN THE TANK**
- **ALL MODELS HAVE DESICCANT AIR DRYING SYSTEMS**
- **ENERGY EFFICIENT MOTORS REDUCE THE NEED FOR COSTLY BOOST TRANSFORMERS**
- **COMPRESSORS ARE AVAILABLE IN SINGLE, DUAL OR TRIPLE HEAD CONFIGURATIONS**
- **TEFLON COATED PISTONS AND SLEEVES INSURE LONG LIFE OF THE COMPRESSOR HEADS**
- **AIR INTAKE FILTERS TRAP DIRT, DUST AND AIR BORN POLLUTANTS**
- **ENERGY EFFICIENT THERMALLY PROTECTED MOTORS**
- **MOISTURE INDICATORS STANDARD ON ALL MODELS**

 **TECH WEST INC.**


Manufacturers of Dental Vacuum
and Air Systems

2625 N. Argyle Ave. • Fresno, CA 93727
(559) 291-1650 • (800) 428-7139 • FAX (559) 348-9677

OILLESS AIR COMPRESSOR

Tech West Ultra Clean Oilless Air Compressor provides 100% oil free air to the dental office. With our high performance drying system we guarantee the clean dry air you expect. The design of the Tech West Oilless Compressor provides smooth vibration free running at 72 DB saving the expense and hassle of a sound cover. They are offered in a variety of configurations; singles, duals or triples to meet any need.

Tech West will provide full plumbing and electrical design upon request.

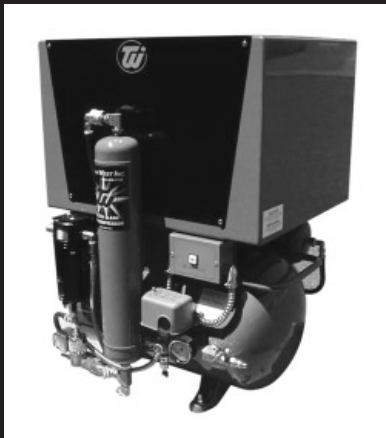

FEATURES

Redesigned Air Dryer makes the unit more compact and reliable.

American made electrical motor with capacitor start and run windings.

Motor will start under full compressor load with a voltage range of 208/248. No Italian or European motor can.

One handle dryer shutdown.

Moisture indicators on the outlet air line.

Integrated carrying handles on the tank.

BENEFITS

This allows the compressor to be installed under cabinets and in closets, saving space.

Faster starting and lower running amps saves electricity.

Compressors can be installed where others would need a buck/boost transformer saving hundreds of dollars.

Quick shutdown and ability to continue working until a service call can be made.

Early indication of moisture in the air before any condensation occurs.

Makes for easy moving and installation.

OILLESS AIR COMPRESSOR

MODEL NUMBER	MAX USERS	AIR DRYER	QUIET COVER	# OF HEADS	TOTAL H.P.	CFM	TANK GALLONS	VOLTAGE	TOTAL AMPS	BREAKER SIZE	DB(A)	EQUIP. WEIGHT	EQUIP. HEIGHT	EQUIP. WIDTH	EQUIP. DEPTH	SHIPPING WEIGHT	SHIPPING HEIGHT	SHIPPING WIDTH	SHIPPING DEPTH
ACO2S1	2	YES	NO	1	1	4.8	10	115	17	20	70	159	27	26	19	189	37	35	30
ACO2S1Q	2	YES	YES	1	1	4.8	10	115	17	20	70	190	31	26	23	220	43	42	31
ACO2S2	2	YES	NO	1	1	4.8	10	230	10	20	70	159	27	26	19	189	37	35	30
ACO2S2Q	2	YES	YES	1	1	4.8	10	230	10	20	70	190	31	26	23	220	43	42	31
ACO4D2	4	YES	NO	2	2	9.6	20	230	20	30	74	228	31	36	21	278	43	42	31
ACO4D2Q	4	YES	YES	2	3	9.6	20	208/230	20	30	70	305	35	37	23	240	43	42	31
ACO6T2	6	YES	NO	3	3	14.6	30	230	30	40	77	298	33	41	24	348	43	42	31

Revised 1/08

**ALL OILLESS COMPRESSORS CARRY A FULL THREE YEAR WARRANTY ON OILLESS HEADS
AND A FULL TWO YEAR WARRANTY ON ALL OTHER COMPONENTS**

Exceeds NFPA 99C Guidelines



C 214014 US

TECH WEST INC.
Manufacturers of Dental Vacuum
and Air Systems

2625 N. Argyle Ave. • Fresno, CA 93727
(559) 291-1650 • (800) 428-7139 • FAX (559) 348-9677

FDA
RECOGNIZED